



Extracorporeal Shock Wave Therapy

# ZEUS WAVE

Radial(Pneumatic-type) Pain Treatment System



Well-being  
With  
Wever



**Wever Instruments Co., Ltd.**

Address Wever B/D, 22-9, Chungui-ro 37beon-gil, Uijeongbu-si Gyeonggi-do, Republic of Korea  
Tel +82-31-853-1924 Fax +82-31-853-1925 Website www.weverinstruments.com  
ZEUS WAVE(Pneumatic-type) PR brochure ENG Rev.1




# ZEUS WAVE

Extracorporeal Shock Wave Therapy

## ○ Pneumatic Type

### What is ESWT

Extracorporeal shock wave therapy (ESWT) is one of the few proven medical techniques that can safely treat, in a non-surgical/non-invasive way, acute and chronic pain syndromes in the musculoskeletal system, on which traditional pain treatment had not shown a great effect.

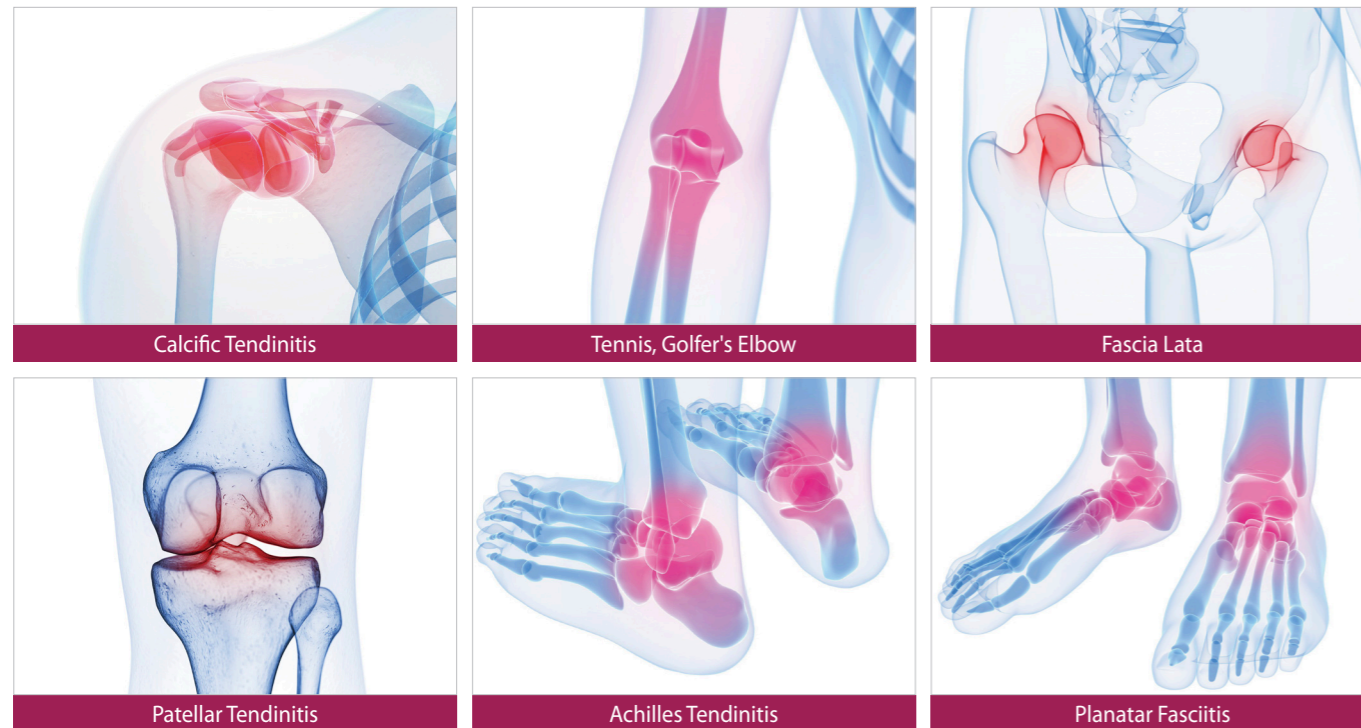


### Application of FSWT

When treating enthesopathies (Tendon, Bone)	When breaking up calcium deposits	When identifying the location of the pain spots with the pain triggers that cause synalgia	When treating superficial pain trigger points / Pain points and profundus pain trigger points

### Principle of Radial

When the lesion site is repeatedly irradiated with shock wave energy generated from the gun head, temporary micro-damage to the lesion site induces the formation of neovascularization, and the supply of blood flow increases to facilitate the healing of the lesion tissue, thereby relieving the pain. This mode is effective in relieving pain by destroying the calcium deposits accumulated at the trigger point.

### Features of RSWT

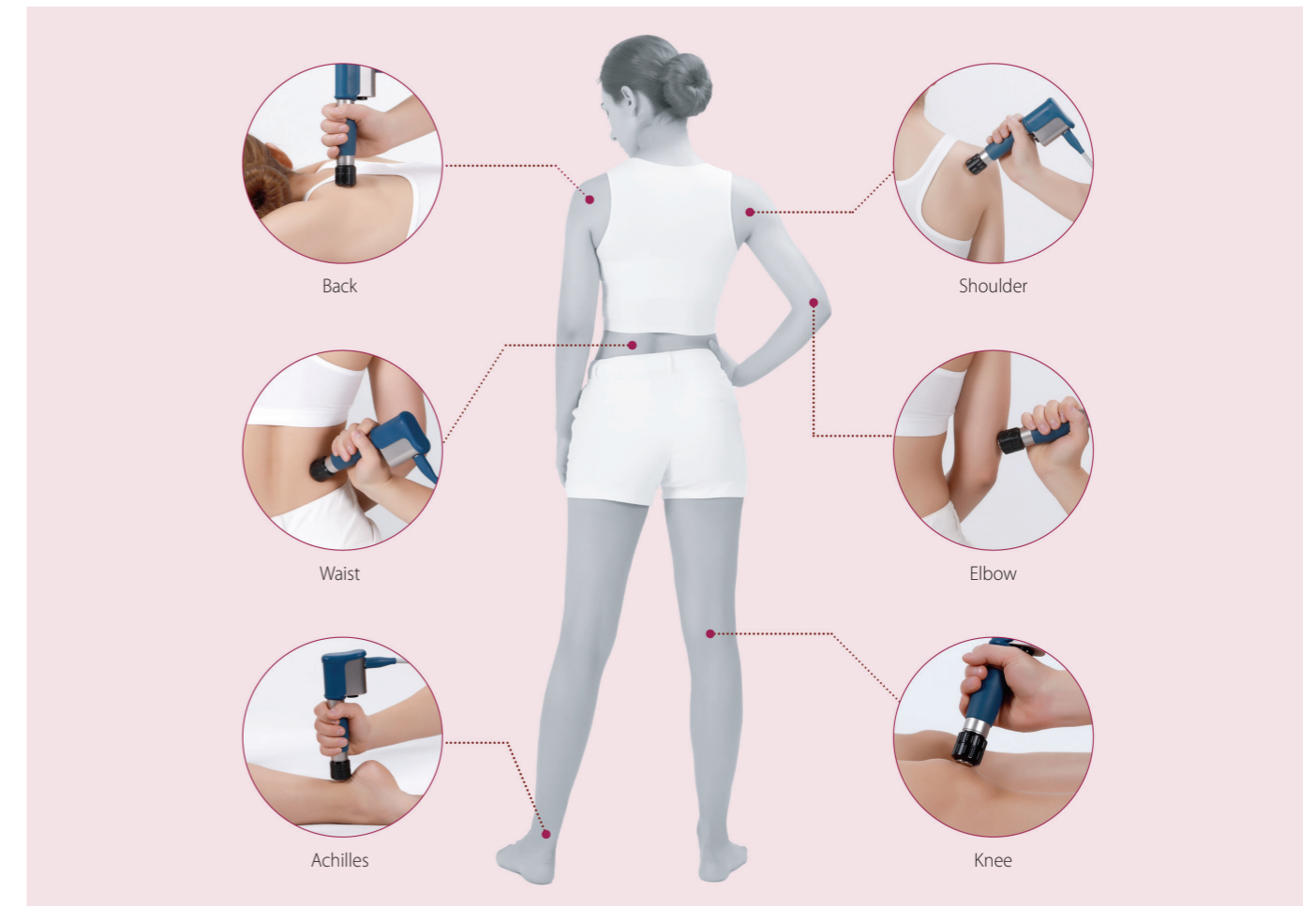

### Air-gun Handpiece & Applicators



### ZEUS WAVE Pneumatic-type Specifications

<b>Rated Power</b>	AC 200~240V	<b>Frequency</b>	0~10Hz ± 10%
<b>Rated Frequency</b>	50/60Hz	<b>No. of shocks</b>	0~500 Shock ± 10%
<b>Power Consumption</b>	1.1kVA	<b>Dimension</b>	500(W) x 720(D) x 980(H)mm
<b>Max. intensity</b>	7.0bar ± 10%	<b>Weight</b>	About 65kg